

red<mark>dot</mark> award 2017 winner

## FENNIAPRIZE 17

GRAND PRIX

# Q-Flow — Award-winning surgical light




# Fluent usability

## Since 1901


"In addition to the exquisite quality of work, we have taken the utmost interest in the functionality of the equipment. In this regard, several doctors committed to the cause have contributed their advice, based on their valuable experience."

Juho Merivaara in 1926, founder of Merivaara

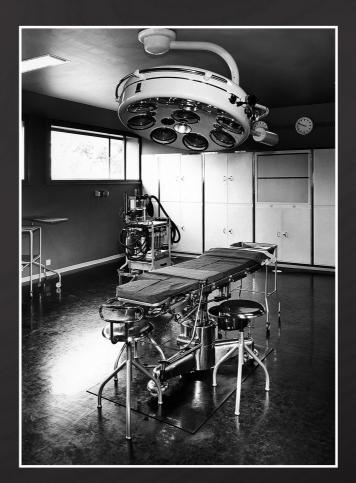

In a modern operating room from the 1950s tables and lights were manually adjusted.




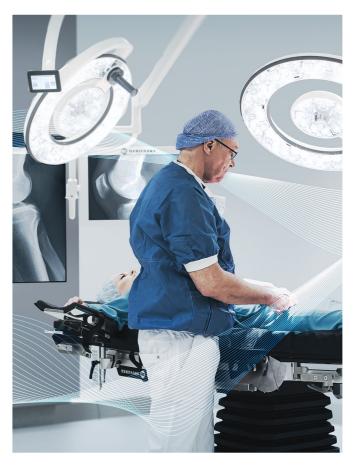

We are proud that the company's values 120 years ago were so ahead of their time that they can still today be fully applied when designing new innovative solutions for healthcare professionals. This is demonstrated by our latest award-winning products, the Q-Flow surgical light and the Smarter Practico operating table with functionality and ease of use at the heart of the design. Our user interfaces are developed together with hospital personnel so that all our products and systems can be used intuitively.

We call this Fluent Usability; operating room teams can focus on patient care, not the management of complex technologies. Top surgeons appreciate this all over the world.



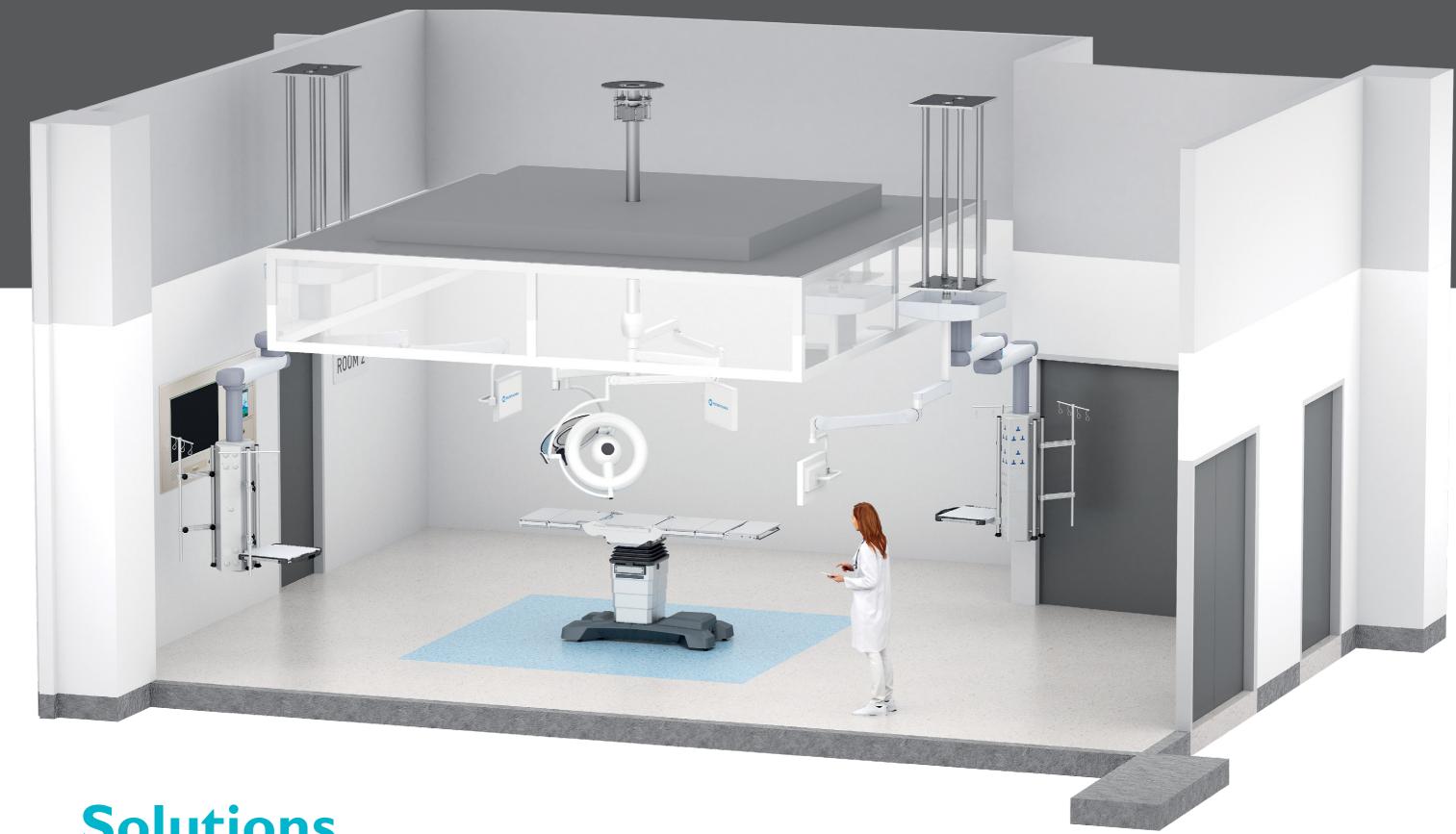

## **Solutions**

Merivaara Solutions offer operating rooms a wide range of high-quality products, systems and services. Our solutions not only guarantee patient safety, but also provide healthcare operators with the best value for money.

When designing our solutions around the entire surgical team, it is important for us to listen to the customer's needs. This allows us to ensure that the team can implement the procedures safely, quickly and reliably. At the same time, workflow is optimised and the productivity of the entire team will increase.

Merivaara has long experience from leading and implementing projects in operating rooms. For us, keeping projects to schedule and reducing costs is extremely important.

• SURGICAL LIGHTS • EXAMINATION LIGHTS • OPERATING TABLES • TABLE ACCESSORIES INTEGRATION SYSTEMS • MONITORS & MONITOR ARMS • PENDANTS • PROJECTS • MERIVAARA SERVICES

## Benefits of the Q-Flow family



The design of the Q-Flow light is optimised to improve air flow circulation in the operating area. This will help keep the area clean, improving working conditions for the entire surgical team and increasing patient safety. Turbulence intensity is remarkably low 15.9%, when the industry standard specifies it to be below 37.5 %.



Merivaara's patented Intueri<sup>™</sup> system provides surgeons with the ability to adjust the dimming and the field size of the light field from the sterile area, allowing surgeons to keep their focus on the operating area. The Intuitive user interfaces of the Q-Flow touch panel and its auxiliary controllers – multifunctional remote control Merimote<sup>™</sup> and integration system OpenOR<sup>™</sup> – offer similar, simplified controls for the light and its camera, in multiple languages. This allows users to switch between control devices smoothly as needed.

The in-light camera unit is completely hidden behind the cover glass, leaving no obstacle in the surgeon's peripheral view. Thanks to the central position of the camera unit, the camera can always be centred in the operating range. Wireless technology enables the suspension arms to rotate a full 360 degrees, allowing surgeons to focus the camera at an optimal viewing angle.

## Design optimised for OR's air flow

Hospital associated infections (HAIs) are a highly recognised challenge in hospitals around the world. HAIs increase morbidity, mortality and length of hospital stay, adding to health care costs. More than 4 million patients in Europe and approximately 1.7 million in the US are affected annually. The prevalence of HAIs in Europe is around 7.1 %.

One cause of contamination is uncontrolled air flow in operating rooms. Q-Flow<sup>™</sup> is designed for an optimised air flow that allows ventilation to work properly in an operating room. Traditional surgical light heads cause air to rise in the operating area, resulting in increased particle content and an increased infection burden for the patient. Due to the optimally designed Q-Flow<sup>™</sup>, with a turbulence intensity of only 15.9 %, there is no additional particle burden created in the operating area. In addition, this also improves the working conditions for the surgical team, as it helps to keep the area clean from harmful smoke and gases.









### I. REDUCED RISK OF INFECTION

#### 2. PERFECT TISSUE COLOUR RENDERING

The R9 value (red) of the Q-Flow light makes it easy for the surgeon to separate between different tissue and vascular colours.

The RI3 value (skin) 98 is also excellent, which is particularly

important in plastic surgery and in operations where seeing the skin accurately is important.

#### 3. COMFORTABLE WORKING CONDITIONS

Merivaara's Dynamic Obstacle Compensation (DOC<sup>™</sup>) decreases the need to adjust the luminaire.

• Shadowless light in the surgical area provides optimal light in all circumstances

• Automatic light dimming reduces temperature above surgical team's

The light restores original lighting conditions automatically

### 4. INTUITIVE TO USE

### 5. ERGONOMIC CAMERA SOLUTIONS



# **Q-Flow models**

### Fluent

Main features:

- Excellent optical performance
- Design optimised for air flow
- Easy to clean
- MeriMote and
- OpenOR compatible

## Vision

Additional features:

- HD camera option
- Green ambilite
- Lightweight


## Flexible viewing angles

## Intelligent


Additional features:

## Maximum comfort and ease of use

## Functional simplicity

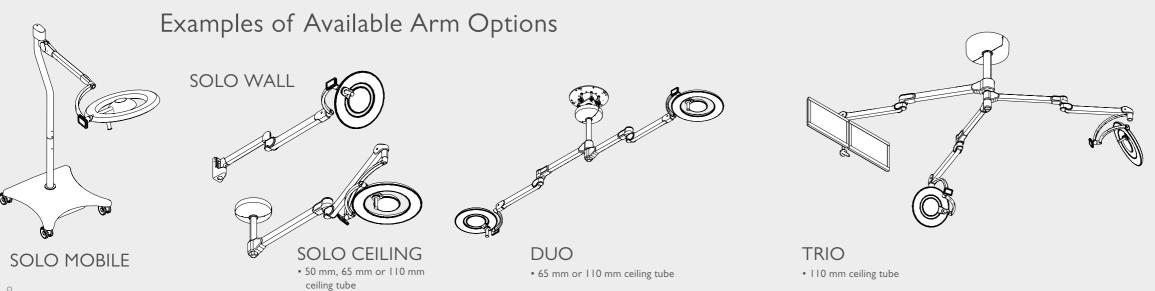

• Intueri<sup>™</sup> Sterile user interface • DOC<sup>™</sup> Dynamic Obstacle compensation

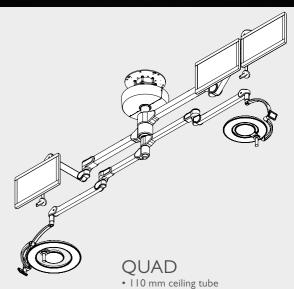

## Our solutions for Q-Flow



#### OpenOR integration system

OpenOR enables you to connect all video and audio sources, medical devices, and room functions to be displayed on monitors or info screens. OpenOR can also control lights, tables, ventilation and blinds, and it communicates with hospital information systems (HIS) and building management systems (BMS).

## Accessories



Disposable sterile light handle cover A43460601 (Set of 200 units)



Display protectors

Provides additional impact protection to your valuable arm mounted medical

clean and maintain. Installed at our factory prior to shipment.

grade monitors. Optically fully transparent surface reduces glare for optimal image

visibility and prevents scratches. Can be removed, cleaned and reapplied. Easy to

Video transfer solutions

to other devices in the hospital.

We can offer complete solutions for transferring live video stream from the Q-Flow lamp camera

ight Sterilisable Handle A43280901 Kit for Q-Flow Lamp Head (10 handles)



MeriMote 5202691 Remote control unit with wall mounting kit Multifunctional remote for Q-Flow systems • 2-way wireless control for all Q-Flow lamp heads • Full in-light camera control • Synchronised control of several lamp heads simultaneously



Wireless HD Camera Set 5202823 Superb high-definition image quality with optimised ergonomics

Fully recessed camera unit with rotation possibility

• Centralised positioning, always focused on to the surgical site

• Integrated camera control panel on the lamp head



Monitors & Monitor arms Merivaara offers customised monitor arms for medical monitor mounting.

## **Technical Specifications**

	FLUENT		VISION		INTELLIGENT	
	6F	4F	6	4	<b>6</b> i	<b>4i</b>
Max illumination @ 1.0 m	160 klux	140 klux	160 klux	140 klux	160 klux	140 klux
Colour rendering index (Ra)*	98	98	98	98	98	98
Red colour rendering index (R9)*	98	98	98	98	98	98
"Skin colour" rendering index (RI3)*	98	98	98	98	98	98
Turbulence intensity, DIN 1946	15.9 %	<35 %	15.9 %	<35 %	15.9 %	<35 %
Variable colour temperatures (K)**	3700 / 4400 / 5100					
Depth of illumination (L1+L2) @ 60 %	1200 mm	690 mm	1200 mm	690 mm	1200 mm	690 mm
Depth of illumination (L1+L2) @ 20 %	1800 mm	1500 mm	1800 mm	1500 mm	1800 mm	1500 mm
Light field diameter	200–370 mm	200–320 mm	200–370 mm	200–320 mm	200–370 mm	200–320 mm
Integrated dimming	10-100 %	10-100 %	10–100 %	10–100 %	10–100 %	10-100 %
Number of LEDs	90	69	90	69	90	69
Average life time of LED	>60 000 h					
External dimensions (diameter)	700 mm	560 mm	700 mm	560 mm	700 mm	560 mm
Lamp weight	16 kg	13 kg	16 kg	13 kg	16 kg	13 kg
Ambilite	White	White	Green	Green	Green	Green
Touch panel control	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$
OpenOR™ compatible	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$
MeriMote <sup>™</sup> compatible	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$
Extended autonomy	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$
Integrated camera control (touch screen)	-	-	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$
Wireless full HD camera (1080p60)	-	-	0	0	0	0
Detachable sterilisable handle	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$
Disposable sterile handle cover	0	0	0	0	0	0
Intuitive Sterile Surgeon Control (Intueri™)	-	-	-	-	$\checkmark$	$\checkmark$
Fluent Sterile Surgeon Control		$\checkmark$	-	-	-	-
Shadowless optical system		$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$	$\checkmark$
Dynamic Obstacle Compensation (DOC™)	-	-	-	-	$\checkmark$	$\checkmark$

 $\sqrt{1}$  = included O = optional

\*) Tolerance ±3% \*\*) Tolerance ±200 °K

0 |

### All values are measured in test laboratory conditions EN $60601\mathchar`-2\mathchar`-41$



### Health Technology with a Human Touch

Made in Finland with respect for our unique nature.





Distributed in the United Kingdom and Ireland by Bender UK Low Mill Business Park Ulverston Cumbria LA12 9EE

Tel: +44(0)1229 480123 Tel ROI: +353 1 5060611

internalsales@bender-uk.com

www.bender-uk.com

The status of the information, specifications and illustrations in this publication is indicated by the date of publication. Whilst the greatest care has been taken over its compilation, Merivaara Corp. is not responsible for any errors or omissions. Merivaara Corp. reserves the right to make changes to the technology, features, specifications and design of the equipment and models without notice. Please check the Merivaara website for the latest updated specifications and other information.

All trademarks are Merivaara property unless otherwise stated. Patents granted or pending.

©Merivaara Corp. DOC ID: 6099EN-9, Date: 7/2022 \_\_\_\_\_